

80*67mm

BIONIME
BIONIME CORPORATION
No.100, Sec. 2, Daqing St., South Dist.,
Taichung City 40242, Taiwan
P. +886 4 2369 2388
F. +886 4 2261 7586
Email: info@bionime.com



RIGHTEST
Official Website

101-3HM510-XXX EN(SG) Rev. Date: 2026-09

RIGHTEST™ | ifree² cgm CONTINUOUS GLUCOSE
MONITORING SYSTEM

USER'S MANUAL

USE WITH RIGHTEST IFREE CGM APP

**READ THIS FIRST:**

It is important to read the entire contents of this manual before using the RIGHTEST iFree 2 Continuous Glucose Monitoring System. The instructions, warnings, precautions, safety information and tips contained within this manual are intended to ensure proper use and optimal results. Discuss the best way to use your RIGHTEST iFree 2 CGM with your healthcare professional. Failure to operate the system according to the guidelines and safeguards specified in this manual may present risks. If your glucose readings do not match your symptoms or how you are feeling, check your fingerstick blood glucose level with a blood glucose meter and consult a healthcare professional if necessary.

Indications For Use & Statement And Advisory	6
Intended use	6
Indications For Use	6
Statement And Advisory	6
Intended Benefit	7
Potential Side Effect	7
Contraindication	7
Safety Information	9
Symbol Definitions	15
Getting To Know Your iFree 2 CGM	18
Performance Features	18
Safety Features	18
CGMs Components	20

Accessories	22
Before You Start	22
Install The Mobile App	22
Setting Up Your iFree 2 CGM	24
Scan The Sensor Kit	24
Apply Your Sensor	26
Attach Your Transmitter	35
Connect Transmitter With Display Device	37
Ending A Monitoring Session	38
End The Monitoring Session	38
Sensor And Transmitter Removal	39
Understand Your Glucose Readings	40
Home Screen Indicator And Displays Overview	40

Glucose Trend Arrow And Arrow Color	42
HIGH/LOW Readings	44
Adding Diaries	44
AGP Report, Insight Report, And Key Indicators	44
Calibration	45
Connection And Data Uploading	46
Treatment Decisions	47
When Not To Use Sensor Readings To Make Treatment Decisions	49
Trend Arrows And Treatment Decisions	50
Specifications	55
Sensor Kit Specifications	55
Transmitter Specifications	57
IT Networks Characteristics And IT Security Measures For App	59

Appendix	60
Alert Audible Output	60
System Alerts	60
Customer Service	71
Device Performance Characteristics	82

Indications For Use & Statement And Advisory

Intended use

The RIGHTEST iFree 2 Continuous Glucose Monitoring System (hereafter referred to as the "iFree 2 CGM") is indicated for detecting glycemic trends and for the management of diabetes in persons aged 18 and older. It is designed to replace fingerstick blood glucose testing for diabetes treatment decisions. Interpretation of the iFree 2 CGM results is based on the glucose trends and several sequential readings over time. It also aids in detecting episodes of hyperglycemia and hypoglycemia, facilitating both acute and long-term therapy adjustments.

Indications For Use

For persons aged 18 and older with type 1 or type 2 diabetes who use insulin or other diabetes medications requiring regular glucose monitoring.

Statement And Advisory

This manual is designed to instruct all personnel responsible on the proper use and care of the iFree 2 CGM in non-professional environments. All users are encouraged to read this manual carefully before using the system.

Intended Benefit

Supports user awareness of glucose levels and trends, enabling informed diabetes self-management.

Potential Side Effect

Skin Irritation Reaction Caused by the Sensor Adhesive: Some individuals may be sensitive to the medical adhesive that keeps the sensor attached to the skin. Skin reactions such as redness, bleeding, inflammation, bruising, or itching may occur at the sensor insertion site. Remove the sensor if these occur while wearing it; if they occur after removal, monitor the site and consult your healthcare professional if symptoms persist or worsen. If necessary, consult your healthcare professional.

Contraindication

- Known allergy to disinfecting alcohol or medical adhesives.
- The iFree 2 CGM is not recommended for insertion into skin areas affected by tattoos, scars, burns, wounds, inflammation, or extensive lesions, including psoriasis, eczema, and dermatitis herpetiformis, as these conditions may interfere with proper device placement and compromise measurement accuracy. The device

should be used with caution in patients with sensitive skin, a predisposition to skin ulceration, known allergies to disinfectants or medical adhesives, or an increased risk of bleeding, as these may affect device placement and measurement accuracy.

-  **No MRI/CT/Diathermy:** The iFree 2 CGM (sensor, transmitter, and/or other display devices) must be removed prior to Magnetic Resonance Imaging (MRI), Computed Tomography (CT), or high frequency electrical heat (diathermy) treatment. People who are unable or unwilling to contact with their healthcare professional are not recommended to use the iFree 2 CGM. Sufficient vision or hearing is critical for successful use of the system including effective recognition of the alerts.
- **DO NOT Use If You Are Pregnant, on Dialysis or Critically Ill.** Do not use the iFree 2 CGM if you are pregnant, on dialysis or critically ill; or on users with other implanted medical devices (e.g., a pacemaker). The system has not been evaluated for use in these populations.
- iFree 2 CGM is contraindicated in patients with infection, edema at or near the sensor insertion site, those receiving vasoactive drug therapy, or with impaired tissue perfusion.
- iFree 2 CGM should not be used during episodes of severe

hypoglycemia (<40 mg/dL), severe hyperglycemia (>500 mg/dL), or during periods of rapid blood glucose fluctuation.

Safety Information

The following is a summary of safety information which must be observed before using the iFree 2 CGM. WARNING indicates potential danger to the user. PRECAUTION indicates potential injury to the user or damage to the system. To minimize risks, read the following safety information before using the system. Improper use and maintenance may damage the system resulting in failure or injury to the user. It is important to understand that this safety information is not exhaustive. It is intended to ensure the safety of the user when using the system.

If you experience serious incidents caused by the use of the iFree 2 CGM, contact your local emergency services for help. Please report the incident to Bionime Corporation and the local competent authority.



WARNINGS:

Use a Blood Glucose (BG) Meter To Make Treatment Decisions Under The Following Conditions:

- During the 1-hour warmup period when you start a new sensor. You will not receive any sensor readings, or alerts until your system begins to transmit data.
- If you suspect that your sensor readings may be inaccurate for any reason.
- If your sensor readings do not match what you are feeling.
- If you are experiencing symptoms that may be due to low or high blood glucose.
- If your system does not include your current glucose concentration or a glucose trend arrow.
- If you wish to confirm hypoglycemia or impending hypoglycemia as reported by the system.
- If you are experiencing rapid glucose changes (greater than 2 mg/dL per minute), the sensor readings displayed may be less accurate and less timely.
- If your system shows "Use Blood Glucose Meter (BGM) .

Not Receiving Urgent Alerts or Alerts under the Following Conditions:

- When either your display device or transmitter battery is dead.

- When your display device is turned off or when RIGHTTEST iFree CGM App is not running in the background.
- When there is a system error (e.g., no glucose readings, sensor error, signal loss, etc.) or damage to the system.
- During the 1-hour sensor warmup period.
- If your glucose readings do not match how you are feeling.
- When the display device is out of range from your transmitter; or obstacles (metal, walls, water, etc.) are between them; or when you are in crowded place (station, exhibition etc.).

Modification of the System is Not Permitted: Do not modify or tamper with any components or accessories of the iFree 2 CGM. Do not use any component of the iFree 2 CGM with any product not included in this system. Otherwise, you may damage the integrity of the system and put yourself at risk especially when you have a severe low or high glucose event.

Do not allow children or pets to play with: Do not allow children or pets to play with any parts of iFree 2 CGM.

Choking Hazard: The iFree 2 CGM contains small components that may be dangerous if swallowed.



PRECAUTIONS:

Calibration Safety: iFree 2 CGM is factory calibrated and ready to use after 1-hour warmup period without additional calibration by user. During monitoring session, calibrating the iFree 2 CGM is optional. If you choose to calibrate, use only fingerstick blood glucose values within 5 minutes to align the CGM readings with them. Entering incorrect values, or using blood glucose readings obtained from alternate sites, may result in inaccurate sensor readings and could lead to missed high or low glucose events.

Skin Irritation Reaction Caused by the Sensor Adhesive: Some individuals may be sensitive to the medical adhesive that keeps the sensor attached to the skin. If you develop a rash around or under your sensor, remove the sensor and stop using the iFree 2 CGM. If necessary, consult your healthcare professional.

Avoid Skin Care Products: Do not apply skin care products such as sunscreen, moisturizer, perfume or insect repellent over the sensor insertion site (less than 1 inch or 2.5 cm) or any components of the iFree 2 CGM. Failure to comply may

lead to damage of the plastic used in the iFree 2 CGM or reduction in the stickiness of the sensor adhesive.

Do NOT Attempt to Reinsert a Sensor: If the adhesive patch is loose or if the sensor tip is pulled out from your skin, remove the sensor and replace it with a new sensor. Sensor readings may be unreliable until a new sensor is inserted.

DO NOT Expose to Extreme Conditions: Avoid freezing the sensor or exposing the components to direct sunlight and environments outside the specified ranges. Exposing these components to out-of-specification conditions may damage the device and result in no or inaccurate glucose readings.

DO NOT Reuse Your Sensor or Inserter: The entire Sensor Kit package is sterilized and designed for single use. It is not suitable for re-sterilization. Re-sterilization of these components may result in no glucose readings and infections.

DO NOT Use If Any Component Appears to be Damaged: A damaged or cracked sensor kit or transmitter may compromise the integrity of the system and contribute to infection risk.

Traveling by Air: Always check and follow flight rules and

regulations before departure. Notify the security personnel of the presence of the iFree 2 CGM and comply with requirements for pat-downs, visual inspection and metal detectors. You must comply with any requests by airline personnel (e.g., turning off the system). Do not pass through an advanced imaging technology (AIT) body scanner (e.g., millimeter wave scanners) or put iFree 2 CGM components through X-ray machines since the effect of this equipment on iFree 2 CGM has not been evaluated.

Changing Date/Time Is Not Permitted: You are not allowed to change date and time during the 15-day monitoring period. Changing the time or date settings during monitoring may result in gaps in the graph or hidden glucose readings.

Keep an Emergency Kit with You: Make sure necessary supplies are always available. Let your family, co-workers, or friends know where the emergency kit is.

The emergency kit should contain:

- Fast-acting glucose tablets.
- Blood glucose monitoring supplies.
- Insulin syringe and rapid-acting insulin (with dosage

instructions from your healthcare professional).

- Adhesive dressing.
- Glucagon emergency kit.

Troubleshooting: If any situation not mentioned in this user manual occurs, please contact your healthcare professional or Customer Service. Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

Be Careful of Electromagnetic Disturbance: Stacking equipment may negatively influence on electromagnetic compatibility. Stay a distance greater than 30 cm (12 inches) from any part of any portable RF communications equipment and at least 1 meter from sensitive equipment. If abnormal behavior is observed due to EM disturbances, please relocate the device accordingly.

Symbol Definitions



Use-by date



Date of manufacture

	Manufacturer		Do not re-use
	Importer		Batch code
	Serial number		Temperature limitation
	Direct current		MR Unsafe
	Caution		Keep away from sunlight
	Humidity limitation		Catalogue number
	Medical device		Unique device identifier
	Keep dry		Sterilized using irradiation
	Type BF applied part device		

	Do not use if the package is damaged
	Discard this product according to local regulations
	Operating Instructions; consult manual for further instructions
IP21	Protected from touch by fingers and objects greater than 12.5 millimeters. Protected from vertical dripping water for 10 minutes.
IP22	Protected from touch by fingers and objects greater than 12.5 millimeters. Protected from dripping water when tilted at 15° for 10 minutes.
IP48	Protected from tools and small wires greater than 1 millimeter. Protected from immersion 10 feet (3.05 meter) for 30 hours.

Getting To Know Your iFree 2 CGM

Performance Features

The iFree 2 CGM provides glucose readings, trends, and levels every minute.

The iFree 2 CGM has user-friendly features and benefits including:

- Ergonomic design of sensor inserter allows users to insert the sensor safely with just one hand.
- Easy-to-read visual glucose values and trends.
- Powerful transmitter storage holds 15 days of glucose readings with zero data loss.
- Visual and audio alerts for hypoglycemia and hyperglycemia.
- Lightweight sensor and transmitter for maximum comfort.

Safety Features

This system detects the glucose concentration within the interstitial fluid and the blood glucose (BG) meter detects the blood glucose concentration. Therefore, the glucose readings from this system may delay from 5 to 15 minutes comparing with blood glucose values.

The iFree 2 CGM offers a number of important safety features when you use it.

These features include:

- Alerts include visual notification, vibrations and sound, depending on your personalized settings.
- When you are out of your target glucose range, display device alerts you.
- Display device warns you if your glucose level falls below to or below 54 mg/dL.
- Urgent alert settings at 54 mg/dL or below cannot be turned off. Some settings (such as alert volume level or temporarily switching to silent mode) can be adjusted according to your needs.
- Display device notifies you when a sensor has failed, expired or when there are system errors.

CGMs Components

The iFree 2 CGM consists of 3 key parts: a sensor kit, transmitter and display device.



Sensor
Inside

SENSOR KIT

Sensor kit contains with a pre-loaded sensor inside the inserter and does not require user assembly. The inserter helps you place the sensor under your skin with ease. The sensor measures your interstitial glucose level.



TRANSMITTER

The transmitter wirelessly sends your glucose data from the sensor to the display device. The transmitter is for single-use only.



DISPLAY DEVICE - APP INSTALLED SMARTPHONE (iOS or Android) and SMARTWATCH (WatchOS or Wear OS)

The display device provides sensor readings and delivers alerts of high and low glucose readings. The system can be used with a smartphone with the mobile App. You can install the iFree CGM App on your smartphone; if you use a smartwatch in conjunction, you can also view the glucose readings and trends provided by the iFree CGM App on your smartwatch.

Accessories



STORAGE VIAL

The vial is used for storage of transmitter to keep them dry.

Before You Start

Install The Mobile App

iFree CGM App can be downloaded from the Google Play Store or App Store. Start by following the on-screen instructions to complete the initial setup if it is your first time using the App. The screens in this manual may look different from your App because of operating systems or updates, please use the App by following the on-screen instructions. Refer to the original user's manual of your smartphone to learn how to change relevant settings.

Before starting monitoring, please confirm the following settings

below:

- **Bluetooth on and location permission agreed:** Connection between your transmitter and the App is via Bluetooth and location; you will receive sensor readings or alerts after enabling and agreeing to the App permission.
- **Notifications on:** Enable and allow notifications to show on your locked screen.
- **Keep the battery charged:** The App will continue working in background and keep draining your battery, make sure to keep the power sufficient (You may regularly remove other App running the the background of your smartphone).
- **Smartphone powered on and running:** Open the App again if you restart your smartphone.
- **Update manually:** Update the operating systems or Applications automatically may change settings or shut down the App. Always update manually and verify the setting afterward.
- **Do not change the time:** You are not allowed to change the date/time during the 15-day monitoring session. It may result in gaps in the graph or hidden glucose readings.

iFree CGM App is only compatible with certain smartphones and operating systems, please check the official website or contact

customer service for more information about compatible devices. During the monitoring session, the following situations may affect data integrity and glucose readings or cause loss of monitoring data:

1. Logging out of iFree CGM App
2. Using the same phone number to log in to iFree CGM App on a different mobile device
3. Changing the date and time of the mobile device
4. Deleting iFree CGM App

Setting Up Your iFree 2 CGM

Before setting up your iFree 2 CGM, make sure you have everything you need:

- Sensor Kit
- Transmitter
- Alcohol Wipes
- Display device: only one display device can be used at a time.
- Blood Glucose (BG) Meter: refer to the Blood Glucose (BG) Meter manual or the notifications of your display device.

Note. When using Blood Glucose (BG) Meter, please refer to the Blood Glucose (BG) Meter manual.

Scan The Sensor Kit



PRECAUTIONS:

Scan Before Monitoring: Every time you start a new monitoring session, **scan the QR Code of the sensor kit with your display device.** Each sensor kit has its unique QR Code which is attached on the sealing paper. Not scanning the QR Code or using incorrect QR Code or mixed using components (sensor kit and transmitter) from different packages will cause system failure or wrong glucose readings.

The following steps describe how to start a monitoring session. If you are unable to start a monitoring session, please contact Customer Service.

1. Open the App on your smartphone.
2. Tap **[Let's Start]** to start a new glucose monitoring session.
3. Scan the QR Code of the sensor kit by using the built-in scanner, once it scanned successfully, a checkmark (✓) will appear on the screen.
4. Make sure you follow the steps in the next two sections ("Apply Your Sensor" and "Attach Your Transmitter"). After installing the sensor and transmitter, click **[Connect]**.

Apply Your Sensor

The sensor is pre-loaded inside the inserter. Before applying the sensor to your skin, familiarize yourself with the information in this section.

Overview of the Sensor Kit



PRECAUTIONS:

The circle indicates where the sensor needle is located during inserting. Do not touch this area against any part of your body where you do not want to insert a sensor.

Choose an insertion site on the back of upper arm where there is an adequate amount of subcutaneous fat.



The following areas are preferable for insertion:

1. Skin that stays flat during normal daily activities (without bending or folding).
2. An area unlikely to be bumped, pushed, or lain on while sleeping.



The following should NOT be selected for insertion:


1. Any area of the body other than the back of upper arm.
2. Skin that is painful to touch, is higher than surrounding skin, is crusting or bleeding.
3. Areas directly over muscle, scars, moles, tattoos, irritation, stretch marks, bones, or lumps.



The following is not recommended for insertion:

1. Sites that are too close (less than 1 inch or 2.5 cm) to an insulin injection site or previous sensor insertion site.
2. Areas constrained by clothing or accessories and areas which experience high amounts of movement during exercise so as to avoid accidental sensor removal due to excessive sweat or body movement.

Sensor base is parallel to your upper arm. Follow these steps to apply the sensor under your skin. Correct application of a sensor ensure fully attachment of adhesive patch on your skin and help the sensor stay under your skin for up to 15 days.

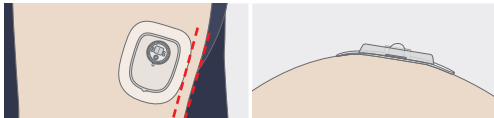

1. Wash and dry your hands.
2. Wipe the insertion site with an alcohol wipe and wait for approximately 2 minutes until the site has fully dried before getting started.

NOTE: Cleaning the insertion site using a plain soap, drying, and then cleaning with an alcohol wipe before insertion of a sensor helps remove any oily residue to let the sensor stick properly. If needed, consider shaving the insertion area to help the sensor stick properly.

3. Open the sensor kit package by peeling off the the sealing paper completely. Take out the Sensor Kit from its package and save the package until the end of the monitoring session.



PRECAUTIONS:

Check the Package. Do not insert the sensor if the sterile package is damaged, broken, or unsealed before you open the package, due to infection risks.


Check the Expiry Date. Discard and do NOT use the Sensor Kit after the expiry date (YYYY-MM-DD) printed on the sealing paper.

4. Peel off the plastic film outside the Sensor Kit.

NOTE: The plastic film prevents the cap opening during transportation.




5. Open the Sensor Kit cap to peel off the release liner from the adhesive patch.

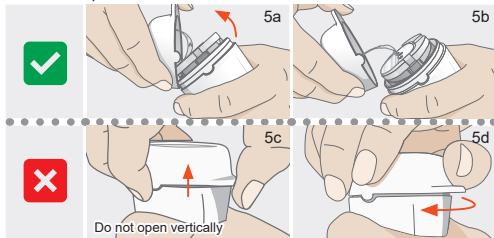

6. Place the inserter over the desired site and push down firmly to insert the sensor. Keep pressing for 10 seconds to ensure the adhesive patch is fully attached to your skin.

NOTE: Rotate the arrow mark toward the shoulder when positioning the inserter to ensure the sensor is secure and comfortable during the wear period. If you are having difficulty inserting the sensor onto the back of your upper arm by yourself,

ask someone to help you or use a mirror for assistance.

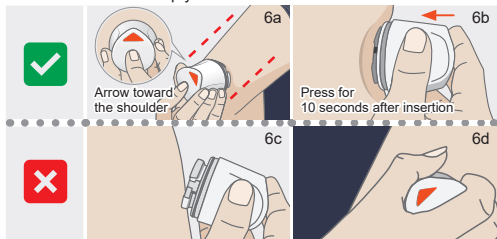


PRECAUTIONS:

1. Apply the sensor immediately after opening its package and the cap. Otherwise, it may present an infection risk.
2. Do not push down the inserter until it is placed over the insertion site.
3. If the insertion is not successful or causes any discomfort, please consult your healthcare professional and use a new

sensor.

4. Do not apply the sensor if it falls out of the inserter when opening the cap.
5. Do not apply the inserter if it is misused or mishandled before insertion.
7. Gently move the inserter away from your insertion site.
NOTE: A temporary skin pressure mark may be visible after removal of the inserter due to the application of pressure during insertion. This is normal and typically resolves within a short period of time.
8. Align both notches on the inserter body and the cap to reconnect them. Discard the inserter in an appropriate puncture-proof or biohazard container according to local regulations for sharps and blood-containing components to prevent cross-contamination and ensure safety.
NOTE: Cap the used sensor inserter immediately after use to avoid needle punching during discarding it or when sensor inserter is mistakenly taken by children.

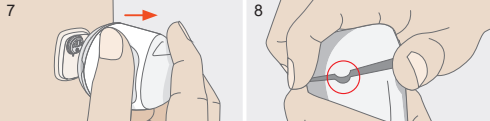


PRECAUTIONS:

Bleeding or bruising at the insertion site under or around the sensor base after applying the sensor is extremely uncommon. If bleeding occurs or you experience high levels of discomfort, follow these steps to reduce risks:

1. Remove the sensor base and adhesive patch. Place sterile gauze or a clean cloth on the bleeding site and apply steady pressure for three minutes or until the bleeding stops.
2. Inspect the site for redness, ongoing bleeding, irritation, pain, tenderness, or inflammation, and contact a healthcare professional for further assistance.
3. Use a new sensor kit and choose another suitable location for reapplication. If the sensor breaks off under the skin, contact a healthcare professional for further assistance.

Overview of the Transmitter

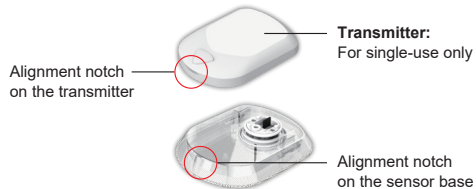

Attach Your Transmitter

Attach your transmitter after the sensor is inserted.

- Make sure the transmitter is stored in the storage vial before use. Do not remove your transmitter until your sensor session is over.



PRECAUTIONS:

Pair Before Use: Make sure the transmitter has been paired with the sensor.

Follow these steps to attach your transmitter:

1. Align the edge of transmitter and the edge of sensor base. Slide transmitter along the edge of sensor base until both notches on the sensor base and transmitter are aligned.

NOTE: After the transmitter and sensor are assembled, they are IP48 rated for water resistance (10 feet or 3.05 meter for 30 hours) and can be worn while bathing, showering, or swimming. Make sure there are no unknown substances on the sensor or sensor base to ensure maximum water resistance.

2. Press down the transmitter until it clicks into the sensor base.

NOTE: Try using a mirror or asking others for assistance to attach your transmitter in the sensor base.

3. Make sure the following instructions:

- (a) All four corners of the transmitter are secured in the sensor base.
- (b) Adhesive patch is fully attached on your skin.



NOTE:

STEP 1: Please ensure that the transmitter is securely attached to the sensor base.

STEP 2: Please "press the pressure-sensitive adhesive patch for 10 seconds daily" to ensure proper adhesion.

Connect Transmitter With Display Device



WARNING:

Use a Blood Glucose (BG) Meter. During the 1-hour sensor warmup period after you insert a new sensor, use a BG meter to make treatment decisions. You will not receive any sensor readings, alerts until your system begins to transmit data.

**WARNING:**

Keep Your Display Device Close. Be sure your display device is close to your transmitter and in the same room. The maximum transmission distance is 6 meters (20 feet) with no obstructions (e.g., walls, metal, glass or water) in between. Obstructions or greater distances may cause Bluetooth signal loss and you may not receive important alerts.

1. Make sure you have followed the steps in the "Apply Your Sensor" and "Attach Your Transmitter" sections. When the screen displays "Confirmation", tap **[Connect]**.

NOTE: Make sure you have installed your sensor and transmitter before you start the following steps.

2. Your smartphone will automatically search for your transmitter. Keep your smartphone close to you.
3. After the system is connected, the screen will display a warmup progress bar. When the warmup is completed, "Warmup" will disappear from your display.

Ending A Monitoring Session**End The Monitoring Session****PRECAUTIONS:**

Do Not Reuse. Reuse of a sensor, sensor base, adhesive patch or transmitter may cause infection or irritation.

Ending a Session Early: If any unexpected issues (irritation or discomfort) happen at the application site, consult your healthcare professional for further assistance.

1. The monitoring session ends automatically when the sensor reaches the ends of its 15-day life and the sensor reading will no longer be shown on the screen. A notification will pop-up to let you know the session has ended. You **MUST** remove or replace the sensor currently in use when you receive this notification. Press **[Remove CGM]** to confirm.
2. To end a monitoring session before receiving the notification, you can select "Stop Monitoring" from "Report". You will see a message warning you that the sensor has not yet expired. Press **[Stop anyway]** to end the session.

Sensor And Transmitter Removal

Do not remove your transmitter until your sensor session is over. Once the session has ended, follow these steps to remove your sensor and transmitter:

1. Grip the edge of the adhesive patch and peel the sensor and transmitter off in one motion.
2. Discard sensor and transmitter according to local regulations for disposal of blood-contacting components.



PRECAUTIONS:

Follow the instructions for proper use of the sensor to avoid sensor breaking.

If the sensor breaks under your skin, do not remove it by yourself. Contact your healthcare professional immediately for further assistance. If any symptoms of infection or inflammation (such as redness, swelling, or pain at the insertion site) occurs, visit medical facility for emergency treatment.

Any remaining adhesive residue on the skin can be removed with warm soapy water or isopropyl alcohol.


Understand Your Glucose Readings

Home Screen Indicator And Displays Overview

Your glucose information (e.g., reading, glucose graph, trend arrow indicating rates and direction of glucose change, etc.) is displayed on the screen of your display device.

It is important to understand these indicators before use.

Overview of Monitoring Screen

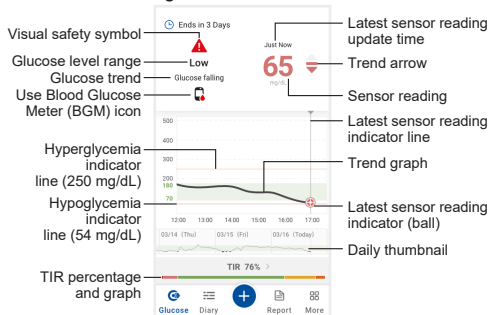

Glucose Trend Arrow And Arrow Color

There are 5 different trend arrows reflecting your glucose readings and how fast they are changing. The color (orange, amber, green, pink and red) of the arrow helps identify the risk of hypoglycaemia and hyperglycaemia.

* "Glucose is steady" means the glucose rate of change is between 0 and 1 mg/dL per minute.

** "Glucose falling/rising" means the glucose rate of change is 1 - 2 mg/dL per minute.

*** "Glucose falling/rising rapidly" means the glucose rate of change is 2 mg/dL per minute or more.

Direction Arrow Color	Glucose is steady*	Glucose rising**	Glucose rising rapidly***	Glucose falling**	Glucose falling rapidly***
Orange > 250 mg/dL					
Amber 181 - 250 mg/dL					
Green 70 - 180 mg/dL					
Pink 54 - 69 mg/dL					
Red < 54 mg/dL					

HIGH/LOW Readings

HIGH means your glucose reading is above 500 mg/dL.

Do a blood glucose test with a BG meter. If you also get a HIGH result (> 500 mg/dL) from the meter test, contact your healthcare practitioner immediately.

LOW means your glucose reading is less than 40 mg/dL. Do a blood glucose test with a BG meter. If you also get a LOW result (< 40 mg/dL) from the meter test, contact your healthcare practitioner immediately.

Adding Diaries

1. Tap the **[+]** icon on the home screen and select **[Add Diary]**.
2. Ensure the diary time is correct.
3. Select categories such as Insulin, Exercise, or Carbs, or utilize the food recognition feature by capturing a photo for automated entry.
4. Follow the on-screen prompts to complete the record, then tap **[Done]** to save.
5. To manage diary records, access the **[Diary]** tab to view, edit, or delete existing entries, or tap a specific diary for detailed information and additional functions.

AGP Report, Insight Report, And Key Indicators

By tapping the **[Report]** icon on the home screen, users can review

the CGM monitoring sessions.

Each monitoring session report includes four key glycemic metrics: Time in Range (TIR), Average Glucose, Glucose Management Indicator (GMI), and Glucose Variability (CV). These metrics allow for a rapid assessment of overall glycemic control and provide detailed information on the frequency of hyperglycemic (high) and hypoglycemic (low) events.

Calibration

The iFree 2 CGM is factory-calibrated and does not require manual calibration. However, manual calibration is optional and can be performed to align the sensor readings with your BG meter values. Refer to the official website to obtain more information.

If you choose to calibrate, take a fingerstick measurement from your BG meter then enter the value according to the following steps:

1. From the Calibrate screen, tap **[Calibrate]**.
2. Enter a exact fingertip blood glucose value measured within 5 minutes then press **[Save]**.
3. You will see a prompt from the screen. Tap **[OK]**.

NOTE:

- Only BG values between 40 mg/dL and 500 mg/dL are

accepted for calibration by the system. If your BG value significantly varied from your sensor reading more than 20%, wait 15 minutes before trying again.

- During periods of rapid glucose change (such as after eating, insulin injection, or exercise), physiological differences between interstitial fluid and blood may cause discrepancies between the sensor reading and the BG meter value. It is advised not to calibrate the CGM during periods of rapid glucose fluctuation to ensure calibration accuracy.

To obtain the best results:

- Calibrate when fasting or before a meal.
- Use only fingertip blood glucose measurements.
- Enter the value immediately after completing the test.
- Do not use previous blood glucose results.
- If calibration fails, wait 15 minutes before trying again.

Connection And Data Uploading

The App can automatically upload your monitoring results to the cloud via the Internet.

Refer to the original manual of your smartphone to learn how to set up a mobile network or Wi-Fi to connect to the Internet. Using a

mobile network to connect to the Internet may be charged for data transmission. It will be charged by your mobile carrier.

Treatment Decisions

Before you start using the iFree 2 CGM for treatment decisions, make sure you are familiar with the tips provided in this chapter and you have a good understanding of how the system works. Refer to the official website to obtain more information.

- Continue to use your blood glucose meter for treatment decisions until you are comfortable with the information you receive.
- Getting familiar with the system could take days, weeks, or even months.
- Work with your healthcare practitioner and follow their recommendations to put together a plan for making treatment decisions including insulin dosing, based on the glucose values provided by the system, to avoid improper timing or incorrect insulin dosage.
- Check your notes to see how carbs, medication, exercise, illness, and stress levels impact your blood glucose readings.



WARNINGS:

DO NOT Ignore Low/High Blood Glucose Symptoms. If your glucose readings do not match how you are feeling, perform a test with a blood glucose meter. Consult your healthcare professional if necessary.

Use a Blood Glucose (BG) Meter To Make Treatment Decisions Under The Following Conditions:

- During the 1-hour warmup period when you start a new sensor. You won't receive any sensor readings, and alerts until your system begins to transmit data.
- If you suspect that your sensor readings may be inaccurate for any reason.
- If your sensor readings do not match what you are feeling.
- If you are experiencing symptoms that may be due to low or high blood glucose.
- If your sensor readings do not include your current glucose concentration or a glucose trend arrow.
- If you wish to confirm hypoglycemia or impending hypoglycemia as reported by the sensor.
- If you are experiencing rapid glucose changes (more than

2 mg/dL per minute), the sensor readings displayed may be less accurate and not as timely.

- If your system shows "Use Blood Glucose Meter (BGM) .

When Not To Use Sensor Readings To Make Treatment Decisions

You must not make treatment decisions based on your sensor glucose reading in the following situations:

- You suspect that your sensor blood glucose readings may be inaccurate for any reason.
- Sensor glucose readings do not match what you are feeling.
- You are experiencing symptoms that may be due to low or high blood glucose.
- The display device shows no glucose information (e.g., an interrupt alert).
- Glucose is Falling/Rising Rapidly (with upward/downward arrow): Glucose measured in interstitial fluid may differ substantially from true blood glucose levels, particularly at times of rapid glucose change (e.g., after meals, taking insulin, or exercising).
- Low Glucose or Urgent Low Message: Sensor glucose readings may not accurately reflect your blood glucose.
- No Glucose Trend Arrow: During the 1-hour warmup period when

you start a new sensor, the system cannot tell you if your glucose is rising quickly or falling quickly.

- **No Current Glucose Concentration and Trend Arrow:** When there is a HIGH/LOW result, you don't have enough information to make a treatment decision.



PRECAUTIONS:

Sensor Readings may be Different from BG Meter Values.

During periods of rapid change in blood glucose (e.g., after eating, taking insulin, or exercising), you may observe differences in glucose readings between interstitial fluid and capillary blood. Due to physiological differences between different body fluids, the sensor readings may be different from fingerstick blood glucose values from BG meters.

Confirm your blood glucose values with a BG meter before making treatment decisions.

Trend Arrows And Treatment Decisions

Trend arrows show the direction and rate of change of your glucose to give you an idea of where your glucose is going. The following table gives you some ideas on how you may use the arrows when considering your treatment.

Trend Arrow	Treatment Decision		
	Low Glucose (< 70 mg/dL)	Glucose in Target Range	High Glucose (> 250 mg/dL)
(No Trend Arrow)	Do a fingerstick blood glucose check with your BG meter. Do NOT treat based on your iFree 2 CGM.		
Double up arrow 	(All Arrow) Do a fingerstick blood glucose check with your BG meter. Do NOT treat based on your iFree 2 CGM.		
Trend Arrow	Treatment Decision		
	Low Glucose (< 70 mg/dL)	Glucose in Target Range	High Glucose (> 250 mg/dL)
Single up arrow 	(Pink/Red Arrow) Do a fingerstick blood glucose check with your BG meter. Do NOT treat based	• If you are about to eat, take insulin to cover your meal. Consider increasing your dose a little since	(Orange Arrow) • If you are about to eat, take insulin to cover your meal. Consider increasing your

	on your iFree 2 CGM.	your glucose is rising. • If you've recently taken insulin or are about to exercise, wait and check your glucose reading later. Avoid "Insulin stacking".	dose a little since your glucose is high and rising. • If you've recently taken insulin or are about to exercise, wait and check your glucose reading later. • If you have not recently taken insulin and have finished exercise, consider adjusting your insulin correction dose upwards. Avoid "Insulin stacking".
--	----------------------	--	--

Trend Arrow	Treatment Decision		
	Low Glucose (< 70 mg/dL)	Glucose in Target Range	High Glucose (> 250 mg/dL)
Horizontal arrow 	(Pink/Red Arrow) Do a fingerstick blood glucose check with your BG meter. Do NOT treat based on your iFree 2 CGM.	• If you are about to eat, take insulin to cover your meal. • If you've recently taken insulin or are about to exercise, wait and check your sensor reading later. Avoid "Insulin stacking".	(Orange Arrow) • If you are about to eat, take insulin to cover your meal. Consider increasing your dose a little since your glucose is high. • If you've recently taken insulin or are about to exercise, wait and check your glucose reading later. • If you have not recently taken

			insulin and have finished exercise, consider adjusting insulin correction dose upwards. Avoid "Insulin stacking".
Trend Arrow	Treatment Decision		
	Low Glucose (< 70 mg/dL)	Glucose in Target Range	High Glucose (> 250 mg/dL)
Single down arrow 	(Pink/Red Arrow) Do a fingerstick blood glucose check with your BG meter. Do NOT treat based on your iFree 2 CGM.	- If you are about to eat, take insulin to cover your meal. Consider taking a lower dose since your glucose is falling. - If you've recently taken insulin or have finished	(Orange Arrow) - If you are about to eat, take insulin to cover your meal. Consider taking a lower dose since your glucose is falling. - If you've recently taken insulin or

		exercise, eat some snacks or fast-acting carbs.	are about to exercise, wait and check your glucose reading later. Avoid “Insulin stacking”.
Trend Arrow	Treatment Decision		
	Low Glucose (< 70 mg/dL)	Glucose in Target Range	High Glucose (> 250 mg/dL)
Double down arrow	 (All Arrow) Do a fingerstick blood glucose check with your BG meter. Do NOT treat based on your iFree 2 CGM.		
Specifications			
Sensor Kit Specifications			
Sensor Glucose Range	40 - 500 mg/dL		
Sensor Use Life	up to 15 days		

Shelf Life	18 months
Sensor Operating Conditions	Temperature: 5°C - 45°C (41°F - 113°F) Relative humidity: 10% - 95%
Sensor Ingress Protection Rating (when installed with a transmitter)	IP48 Protected from tools and small wires greater than 1 millimeter. Protected from immersion 10 feet (3.05 meter) for 30 hours
Storage & Transportation Conditions	Temperature: 5°C - 30°C (41°F - 86°F) Relative humidity: 10% - 90% (in a cool, dry place)
Operating and Storage Altitude	0 to 3,048 meters (0 to 10,000 ft)
Insert Size	52.0 x 60.0 x 61.1 mm (± 0.5 mm)
Sterilization	Sterilized by radiation

Usage	Single use (disposable)
Transmitter Specifications	
Transmitter Size	24.96 x 18.14 x 4.4 mm (± 0.3 mm)
Transmitter Weight	2.4 g with battery (± 0.5g)
Power Source	Single use lithium battery (3.0V, 35 mAh)
Operating Conditions	Temperature: 5°C - 45°C (41°F - 113°F) (During pairing and wearing period.) Relative humidity: 10% - 95%
Operating and Storage Altitude	0 to 3,048 meters (0 to 10,000 ft)
Shelf Life	18 months
Protection Against Electrical Shock	Type BF applied part

Battery Run Time	Up to 15 days
Usage	Single use (disposable)
Storage & Transportation Conditions	Temperature: 5°C - 30°C (41°F - 86°F) Relative humidity: 10% - 90%
Memory Storage	15 days of glucose data (glucose readings stored every minute)
Ingress Protection Rating (when attached to Sensor)	IP48 Protected from tools and small wires greater than 1 millimeter. Protected from immersion 10 feet (3.05 meter) for 30 hours
Data Communication	Bluetooth 4.2 Frequency range BLE: 2402 - 2480 MHz Maximum RF output power of the product: 0 dBm

Data Communication Range	Up to 6 meters (20 feet)
--------------------------	--------------------------

Quality of Service

The transmitter and display device connect to each other via BLE network. Connection quality is in accordance with the Bluetooth Specification v4.2. The iFree 2 CGM System is designed to accept radio frequency (RF) communications from recognized and paired display devices only.

IT Networks Characteristics And IT Security Measures For App

iFree CGM is designed to transmit data between the transmitter and designated display devices.

iFree CGM uses the following interfaces and communication protocols:

Display Device: Bluetooth Low Energy to transmitter. TLS to data platform using cellular data or Wi-Fi. Display Device is only compatible with certain mobile devices and operating system. For using iFree CGM App, please check <https://www.rightest.com/guides> for more information about device compatibility before using

the App.

Use of the iFree CGM requires user registration, and the user should follow instructions on the continuous glucose monitoring system.

Don't pair your CGM over Bluetooth in public areas.

Bluetooth pairing should be done in a private and safe location to reduce cyber risks such as eavesdropping.


In addition to the security provided by the Bluetooth Low Energy connection, communication between the transmitter and mobile applications is protected by additional levels of security and safety mitigations using an encrypted and proprietary data format. This format embeds various industry standard encryption protocols and methods to protect data, verify data integrity, and to detect and prevent data tampering.

Appendix

Alert Audible Output

Based on testing, the maximum alert audible output measured at 100 cm (3 feet) was 62.2 dB(A) when paired with the tested Android device, and 66.4 dB(A) when paired with the tested iOS device (applicable to high and medium priority alerts).

System Alerts

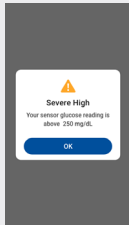

Event - Hyper Alert

Happens when glucose is higher than the set target range.

Visual & Sound Settings - Default Setting:

- Visual Safety Symbol () on Home Screen: Yes.
- Sound & Vibration repeat 3 times or glucose is detected to fall.
- High glucose alert will repeat every 30 minutes or until glucose level returns to the target range.

Editable Setting:

- Sound & Vibration can be turned ON/OFF by the user.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

Event - Hypo Alert

Happens when glucose is lower than set target range.



Visual & Sound Settings - Default Setting:

- Visual Safety Symbol () on Home Screen: Yes.
- Sound & Vibration Alert: Yes.
- Sound & Vibration repeat 3 times or glucose is detected to rise.
- Low glucose alert will repeat every 30 minutes or until glucose level returns to the target range.

Editable Setting:

- Sound & Vibration can be turned ON/OFF by the user.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

Event - Urgent Low Alert

Happens when glucose is at or below 54 mg/dL.




Visual & Sound Settings - Default Setting:

- Visual Safety Symbol () on Home Screen: Yes.
- Sound & Vibration Alert: Yes.
- Sound & Vibration repeat 3 times or until glucose is higher than 54 mg/dL.
- Urgent low glucose alert will repeat every 30 minutes or until glucose level returns to the target range.

Editable Setting:

- Sound & Vibration can be turned ON/OFF by the user.
- Sound volume can be tuned from 10 to 100.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

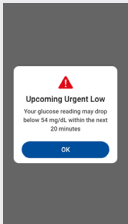

Event - Upcoming Urgent Low

Happens when glucose may drop below the set target value within the next 20 minutes.


Visual & Sound Settings - Default Setting:

- Visual Safety Symbol () on Home Screen: Yes.
- Sound & Vibration Alert: Yes.
- Sound & Vibration repeat 3 times or until the system determines that the readings will not drop below the set target range within the next 20 minutes.
- Upcoming urgent low alert will repeat every 30 minutes or until the system determines that the readings will not drop below the set target value within the next 20 minutes.

Editable Setting:

- Sound & Vibration can be turned ON/OFF by the user.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

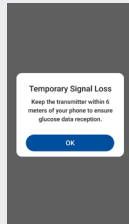

Event - Signal Loss

Transmitter is too far from the display device or when there is an obstacle (e.g., water, wall) in between the transmitter and display device.

Visual & Sound Settings - Default Setting:

- Visual message on the screen: Yes
- Sound & Vibration Alert: Yes.
- Displayed continuously until **[OK]** is pressed.
- System attempts to reconnect every 5 minutes even if user doesn't press **[OK]**
- System alerts will repeat every 30 minutes.


Editable Setting:

- Sound & Vibration can be turned ON/OFF by the user.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

Event - Transmitter Error

Transmitter is not working.

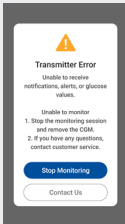


Visual & Sound Settings - Default Setting:

- Visual message on the screen: Yes
- Sound & Vibration Alert: Yes.
- Displayed continuously until **[Stop Monitoring]** is pressed.
- System alerts will repeat every 30 minutes.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

Event - Sensor Fail

The system detects a current error measured by the sensor.

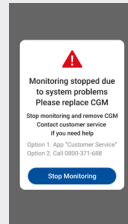


Visual & Sound Settings - Default Setting:

- Visual message on the screen: Yes.
- Sound & Vibration Alert: Yes.
- Displayed continuously until **[Stop Monitoring]** is pressed.
- System alerts will repeat every 30 minutes.

Note: Check your fingerstick blood glucose level with a blood glucose meter.

CGM Data Unstable

The system cannot provide glucose readings temporarily.

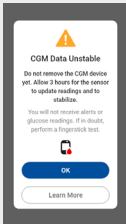


Visual & Sound Settings - Default Setting:

- Visual message on the screen: Yes.
 - Sound & Vibration Alert: No.
 - Displayed continuously until **[OK]** is pressed.
- Note:** Check your fingerstick blood glucose level with a blood glucose meter.

Make Sure Transmitter is Well Attached


Make sure the Transmitter is tightly secured to the base and clicks into place.

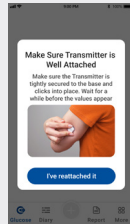

Visual & Sound Settings - Default Setting:

- Visual message on the screen: Yes.
 - Sound & Vibration Alert: No.
 - Displayed continuously until **[I've reattached it]** is pressed.
 - System alerts will repeat every 30 minutes.
- Note:** Check your fingerstick blood glucose level with a blood glucose meter.

Event - Session is Completed
Your sensor session has expired.

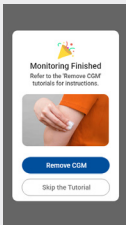

Visual & Sound Settings - Default Setting:

- Visual message on the screen: Yes.
- Sound & Vibration Alert: No.
- Displayed continuously until **[Remove CGM]** is pressed.

Customer Service

We aim to provide great service to our customers. Please review these instructions to make sure you know how to use your product correctly. If you have any questions or encounter any issues with your product, please contact Bionime Customer Service or your authorized distributor.

User's Manual is also available electronically.

- App: "More" tab > Certification > Go to eIFU
- <https://www.rightest.com/guides>
- Free printed copy: Order at <https://www.rightest.com> or contact us.

Tel: +886 4 2369 2388

Fax: +886 4 2261 7586

Email: info@bionime.com

Manufacturer's declaration-electromagnetic emissions

The Transmitter is intended for self-use by the user in the electromagnetic environment (for home healthcare) specified below. The customer or the user of the Transmitter should assure that it is used in such an environment.

Emission test	Compliance	Electromagnetic environment-guidance (for home healthcare environment)
RF emissions CISPR 11	Group 1	The Transmitter uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The Transmitter is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Compliance	

Manufacturer's declaration-electromagnetic immunity			
The Transmitter is intended for self-use by the user in the electromagnetic environment (for home healthcare) specified below. The customer or the user of the Transmitter should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environmentguidance (for home healthcare environment)
Electrostatic discharge (ESD) IEC 61000-4-2	Contact: ± 8 kV Air: ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV	Contact: ± 8 kV Air: ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient/burst IEC 61000-4-4	+ 2 kV for power supply lines, + 1 kV for input/output lines	+ 2 kV for power supply lines Not applicable	Mains power quality should be that of a typical home healthcare environment.

Surge IEC 61000-4-5	+ 0.5 kV, + 1 kV line(s) to line(s), + 0.5 kV, + 1 kV, + 2 kV line(s) to earth	+ 0.5 kV, + 1 kV line(s) to line(s) Not applicable	Mains power quality should be that of a typical home healthcare environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	Voltage dips: 0% UT; 0,5 cycle 0% UT; 1 cycle 70% UT; 25/30 cycles Voltage interruptions: 0% UT; 250/300 cycle	Voltage dips: 0% UT; 0,5 cycle 0% UT; 1 cycle 70% UT; 25 cycles Voltage interruptions: 0% UT; 250 cycles	Mains power quality should be that of a typical home healthcare environment. If the user of the Transmitter requires continued operation during power mains interruptions, it is recommended that the Transmitter be powered from an uninterruptible power supply or a battery.

Power frequency (50, 60 Hz) magnetic field IEC 61000-4-8	30 A/m 50 Hz or 60 Hz	30 A/m 50 Hz	Transmitter power frequency magnetic fields should be at levels characteristic of a typical location in a typical home healthcare environment.
---	--------------------------	-----------------	--

NOTE: UT is the a.c. mains voltage prior to application of the test level

Manufacturer's declaration-electromagnetic immunity			
The Transmitter is intended for self-use by the user in the electromagnetic environment (for home healthcare) specified below. The customer or the user of the Transmitter should assure that it is used in such and environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environmentguidance (for home healthcare environment)
Conducted RF IEC 61000-4-6	3 Vrms: 0,15 MHz - 80 MHz 6 Vrms: in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz	3 Vrms: 0,15 MHz - 80 MHz 6 Vrms: in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz	Portable and mobile RF communications equipment should be used no closer to any part of the Transmitter than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.

Radiated RF IEC 61000-4-3	10 V/m 80 MHz - 2,7 GHz 80% AM at 1 kHz	10 V/m 80 MHz - 2,7 GHz 80% AM at 1 kHz	Recommended separation distance: $d = 1,2\sqrt{P}$ $d = 1,2\sqrt{P}$ 80MHz to 800 MHz $d = 2,3\sqrt{P}$ 800MHz to 2,7 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Interference may occur in the vicinity of equipment marked with the following symbol: 
NOTE1: At 80 MHz and 800 MHz, the higher frequency range applies. NOTE2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

Recommended separation distance between portable and mobile RF communications equipment and the Transmitter

The Transmitter are intended for use in an electromagnetic environment (for home healthcare) in which radiated RF disturbances are controlled. The customer or the user of the Transmitter can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Transmitter as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter (W)	Separation distance according to frequency of transmitter (m)		
	150 kHz to 80 MHz $d = 1,2 \sqrt{P}$	80 MHz to 800 MHz $d = 1,2 \sqrt{P}$	800 MHz to 2,7 GHz $d = 2,3 \sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Manufacturer's declaration-electromagnetic immunity
Test specifications for ENCLOSURE PORT IMMUNITY to RF
wireless communications equipment

The Transmitter is intended for use in the electromagnetic environment (for home healthcare) specified below. The customer or the user of the Transmitter should assure that it is used in such an environment.

Test frequency (MHz)	Band ^{a)} (MHz)	Service ^{a)}	Modulation ^{b)}	Maximum power (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)	Compliance LEVEL (V/m) (for home healthcare)
2450	2400 - 2570	Bluetooth, WLAN, 802.11 b/g/n	Pulse modulation ^{b)} 217 Hz	2	0.3	28	28

NOTE: If necessary to achieve the IMMUNITY TEST LEVEL, the distance between the transmitting antenna and the ME EQUIPMENT or ME SYSTEM may be reduced to 1 m. The 1 m test distance is permitted by IEC 61000-4-3.

- a) For some services, only the uplink frequencies are included.
 b) The carrier shall be modulated using a 50% duty cycle square wave signal.

Manufacturer's declaration-electromagnetic immunity
Test specifications for
ENCLOSURE PORT IMMUNITY to proximity magnetic fields

The Transmitter is intended for use in the electromagnetic environment (for home healthcare) specified below. The customer or the user of the Transmitter should assure that it is used in such an environment.

Frequencies	Test Level [A/m]	Modulation	Dwell time [s]	Compliance LEVEL [A/m] (for home healthcare)
13.56 MHz	7.5	Pulse modulation ^(a) 50 kHz	3	7.5 ^(b)

NOTE:

- (a) The carrier shall be modulated using a 50% duty cycle square wave signal.
 (b) r.m.s., before modulation is applied.

Device Performance Characteristics

The performance characteristics of the iFree 2 CGM system, including accuracy, were evaluated through two clinical studies to assess consistency over time and agreement with the laboratory reference method (YSI 2300). The system demonstrated a MARD within 8.4 - 8.8% under clinical use conditions when compared with the laboratory reference method.

Performance data are shown in mg/dL based on the original test reports and are for technical reference only, not for treatment decisions.

Regression analysis of the Sensors vs. YSI reference Study A

Slope	0.936
Intercept	+11.00 mg/dL
Correlation	0.913
N	9736
Range	40 - 500 mg/dL
Overall mean bias	-1.67 md/dL
Mean Absolute Relative Difference (MARD)	8.8%

Study B

Slope	0.955
Intercept	+3.66 mg/dL
Correlation	0.946
N	5499
Range	40 – 500 mg/dL
Overall mean bias	+2.19 md/dL
Mean Absolute Relative Difference (MARD)	8.4%

Sensor accuracy for all results vs. YSI reference
Study A

Sensor accuracy results for glucose concentrations	Within ± 15 mg/dL (within ± 0.8 mmol/L)	Within ± 20 mg/dL (within ± 1.1 mmol/L)	Within ± 40 mg/dL (within ± 2.2 mmol/L)
< 70 mg/dL	263/282 (93.3%)	274/282 (97.2%)	282/282 (100.0%)
70 - 180 mg/dL	Within $\pm 15\%$	Within $\pm 20\%$	Within $\pm 40\%$
	4837/5936 (81.5%)	5407/5936 (91.1%)	5894/5936 (99.3%)
> 180 mg/dL	Within $\pm 15\%$	Within $\pm 20\%$	Within $\pm 40\%$
	4944/5655 (87.4%)	5310/5655 (93.9%)	5629/5655 (99.5%)
Sensor accuracy results for all results	Within ± 20 mg/dL (within ± 1.1 mmol/L) and within $\pm 20\%$ of reference		
	10965/11873 (92.4%)		

Study B

Sensor accuracy results for glucose concentrations	Within ± 15 mg/dL (within ± 0.8 mmol/L)	Within ± 20 mg/dL (within ± 1.1 mmol/L)	Within ± 40 mg/dL (within ± 2.2 mmol/L)
< 70 mg/dL	361/377 (95.8%)	361/377 (95.8%)	377/377 (100.0%)
70 - 180 mg/dL	Within $\pm 15\%$	Within $\pm 20\%$	Within $\pm 40\%$
	2891/3525 (82.0%)	3220/3525 (91.3%)	3521/3525 (99.9%)
> 180 mg/dL	Within $\pm 15\%$	Within $\pm 20\%$	Within $\pm 40\%$
	1391/1597 (87.1%)	1486/1597 (93.0%)	1597/1597 (100.0%)
Sensor accuracy results for all results	Within ± 20 mg/dL (within ± 1.1 mmol/L) and within $\pm 20\%$ of reference		
	5033/5499 (91.5%)		

Sensor accuracy over wear duration vs. YSI reference Study A

	Beginning (Day 1)	Early Middle (Day 5)	Late Middle (Day 10)	End (Day 14)
Within ± 20 mg/dL (within ± 1.1 mmol/L) and within $\pm 20\%$ of reference	92.9%	95.0%	90.5%	92.2%
Mean Absolute Relative Difference (%)	8.5%	8.1%	9.5%	9.3%

Study B

	Beginning (Day 1)	Early Middle (Day 5)	Late Middle (Day 10)	End (Day 14)
Within ± 20 mg/dL (within ± 1.1 mmol/L) and within $\pm 20\%$ of reference	95.6%	95.6%	86.2%	90.7%
Mean Absolute Relative Difference (%)	7.6%	7.3%	10.1%	8.9%